
beurer

SC 50 Green Planet

DE	Elektrische Zahnbürste Gebrauchsanweisung	6
EN	Electric toothbrush Instructions for use.....	14
FR	Brosse à dents électrique Mode d'emploi.....	21
ES	Cepillo de dientes eléctrico Manual de instrucciones	28
IT	Spazzolino elettrico Istruzioni per l'uso	35
TR	Elektrikli diş fırçası Kullanım kılavuzu.....	42
RU	Электрическая зубная щетка Инструкция по применению	49
PL	Elektryczna szczoteczka do zębów Instrukcja obsługi	56
NL	Elektrische tandenborstel Gebruiksaanwijzing	63
DA	El-tandbørste Betjeningsvejledning.....	70
SV	El-tandborste Bruksanvisning.....	77
NO	Elektrisk tannbørste Bruksanvisning.....	84
FI	Sähköhammasharja Käyttöohje.....	91


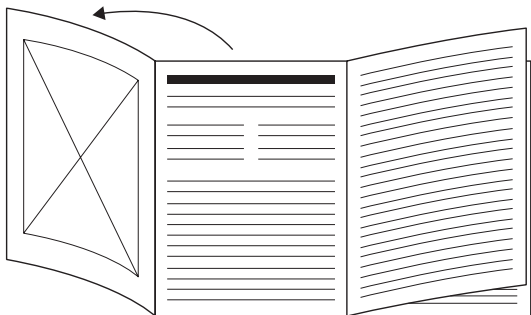
- DE** Klappen Sie vor dem Lesen der Gebrauchsanweisung die Seite 3 aus.
- EN** Unfold page 3 before reading the instructions for use.
- FR** Dépliez la page 3 avant de lire le mode d'emploi.
- ES** Despliegue la página 3 antes de leer las instrucciones de uso.
- IT** Prima di leggere le istruzioni per l'uso aprire la pagina 3.
- TR** Kullanım kılavuzunu okumadan önce 3. sayfayı açın.
- RU** Перед чтением инструкции по применению разложите страницу 3.
- PL** Przed przeczytaniem instrukcji obsługi otworzyć stronę 3.
- NL** Vouw pagina 3 uit om de gebruiksaanwijzing te kunnen lezen.
- DA** Fold side 3 ud, før du læser betjeningsvejledningen.
- SV** Vik ut sidan 3 innan du läser bruksanvisningen.
- NO** Åpne side 3 før du leser bruksanvisningen.
- FI** Käännä sivu 3 auki ennen käyttöohjeen lukemista.


A


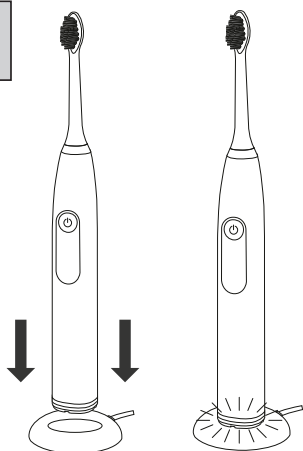
B

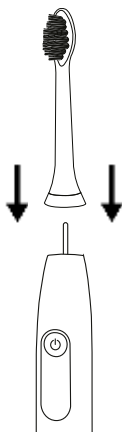
C

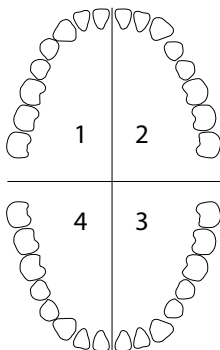D

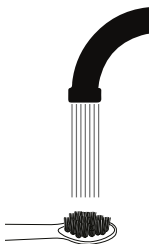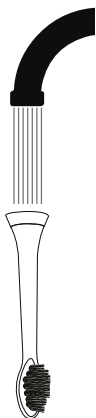E

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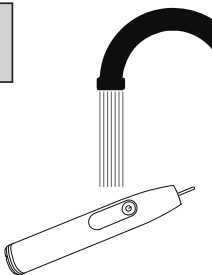Read these instructions for use carefully. Observe the warnings and safety notes. Keep these instructions for use for future reference. Make the instructions for use accessible to other users. If the device is passed on, provide the next user with these instructions for use.

CONTENTS

1. Warnings and safety notes.....	14
2. Included in delivery	16
3. Signs and symbols	16
4. Intended purpose	18
5. Device description.....	18
6. Usage	18
6.1 Preparing for use	18
6.2 Starting use	18
6.3 Travel lock.....	19
7. Cleaning and maintenance.....	19
8. Accessories and replacement parts.....	19
9. Disposal.....	20
10. Technical specifications	20
11. Guarantee	20

1. WARNINGS AND SAFETY NOTES

⚠ WARNING

- Keep the charger and mains adapter (not included in delivery) away from water or other liquids.
- Do not place the charger or mains adapter (not included in delivery) near a sink or bathtub full of water.

- This device may be used by children over the age of 8 and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely, and are fully aware of the consequent risks of use.
- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children unless supervised.
- Do not use the device if it shows signs of damage or does not function properly. In these cases, contact Customer Services.
- If the mains connection cable is damaged, you must dispose of the charger.
- Only use the charger supplied with your device.
- The device may only be charged with a mains adapter  (not included in delivery) with a safety extra-low voltage output of 5 V and max. 1 A.
- Do not pull on the mains cable or on the device to disconnect the mains adapter plug from the socket.
- Under no circumstances should you open or repair the device yourself, as faultless functionality can no longer be guaranteed thereafter. Failure to comply will invalidate the guarantee.

Notes on handling batteries

⚠ WARNING

- If fluid from a battery comes into contact with your skin or eyes, wash the affected areas with water and seek medical assistance.
- Risk of explosion! Do not throw batteries into a fire.
- If a battery has leaked, put on protective gloves and clean the battery compartment with a dry cloth.
- Do not disassemble, split or crush the batteries.

⚠ CAUTION

- Protect batteries from excessive heat.
- Only use chargers specified in the instructions for use.
- Batteries must be charged correctly prior to use. The instructions from the manufacturer and the specifications in these instructions for use regarding correct charging must be observed at all times.
- Fully charge the battery prior to initial use (see Chapter “Preparing for use”).
- In order to achieve as long a battery service life as possible, fully charge the battery at least 2 times a year.

2. INCLUDED IN DELIVERY

Check that the exterior of the cardboard delivery packaging is intact and make sure that all contents are present. Before use, ensure that there is no visible damage to the device or accessories and that all packaging material has been removed. If you have any doubts, do not use the device and contact your retailer or the specified Customer Service address.

- 1 electric toothbrush
- 2 brush heads
- 1 charger with USB cable
- Instructions for use

3. SIGNS AND SYMBOLS

The following symbols are used on the device, in these instructions for use, on the packaging and on the type plate for the device:

⚠ DANGER Indicates an imminent danger. If it is not avoided, death or serious injury will occur.
⚠ WARNING Indicates a potentially imminent hazard. If it is not avoided, there is a risk of death or serious injury.
⚠ CAUTION Indicates a potentially imminent hazard. If it is not avoided, slight or minor injuries may result.
NOTICE Indicates a potentially harmful situation. If it is not avoided, the device or something in its vicinity may be damaged.
 Product information Note on important information

	Read the instructions
	Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment).
	Do not dispose of batteries containing hazardous substances with household waste
	The products demonstrably meet the requirements of the Technical Regulations of the EAEU.
	Standard for independent certification of recycled material
	Packaging made from recycled material
	CE labelling This product satisfies the requirements of the applicable European and national directives.
	Separate the product and packaging elements and dispose of them in accordance with local regulations.
	Marking to identify the packaging material. A = Material code, B = Material number: 1–7 = Plastics, 20–22 = Paper and cardboard
	Conformity assessed marking for the United Kingdom
	Protection class III device
	Importer symbol
	Manufacturer
	Protection against temporary submersion. Device can be immersed to a depth of 1 metre for up to 30 minutes.
 PIN1+ PIN4- 	USB connection
	Direct current

4. INTENDED PURPOSE

Only use the device to clean teeth on humans. Do not use the device on animals! The device is not suitable for children under 8 years. The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage resulting from improper or careless use.

5. DEVICE DESCRIPTION

The associated drawings are shown on page 3.

1 Brush head

2 ON/OFF/Mode button

3 LC display

- 5 modes: Clean, Whiten, Sensitive, Polish, Gum Care
- Emoticons
- Battery level
- Travel lock

4 Light ring

Pressure control

Red light ring: Pressure too high

White light ring: Pressure OK

5 Charger

6. USAGE

6.1 Preparing for use

Charging the electric toothbrush **A**

Fully charge the electric toothbrush before first use. To do so, proceed as follows:

1. Plug the USB charging cable of the inductive charger **5** into a USB port or into a suitable mains adapter with a USB connection  of protection class 2  (not included in delivery) and plug it into a socket.
2. Place the electric toothbrush on the inductive charger **5**. The battery level indicator appears on the LC display **3** whilst the device is charging.

Fitting the brush head **B**

Place the brush head **1** onto the device from above until it clicks into place. Make sure that the brush head **1** points in the same direction as the ON/OFF/Mode button **2**.

Indicator bristles

Indicator bristles lose their colour over time depending on how each person uses the device, thus indicating the best time to change the brush head **1**.

We recommend changing the brush head **1** after 3 months. If the bristles of the brush head are worn or frayed, they may need to be replaced more frequently.

6.2 Starting use

1. To switch on the electric toothbrush, briefly press the ON/OFF/Mode button **2**. The electric toothbrush pauses briefly every 30 seconds to remind you to clean all 4 jaw quadrants for the same amount of time **C**. The electric toothbrush has an automatic timer and vibrates two times after 2 minutes. The toothbrush switches off automatically after 6 minutes.

The electric toothbrush has a pressure sensor. As soon as you press the brush head **1** too firmly onto your teeth or gums during use, the light ring **4** will start to illuminate red.

2. Select your desired mode by briefly pressing the ON/OFF/Mode button **2**. The currently selected mode is shown in the LC display **3** accordingly. The following modes are available:

Daily Clean	For daily cleaning
Whiten	For whiter teeth
Sensitive	For gentle cleaning
Polish	For removing superficial stains
Gum Care	For looking after the gums

3. To cancel the application prematurely, briefly press the ON/OFF/Mode button **2**. If the optimum cleaning time of 2 minutes is not reached, the LC display **3** will show a sad emoticon. Briefly press the ON/OFF/Mode button again to continue cleaning your teeth in the previously selected mode. The timer now continues running from the point that it was stopped.

6.3 Travel lock

The travel lock prevents the device from being switched on by accident (e.g. on trips, when being carried in a suitcase).

1. To activate the travel lock, press and hold the ON/OFF/Mode button **2** for 5 seconds. The travel lock symbol lights up. The travel lock symbol flashes during activation. As soon as the travel lock is activated, the travel lock symbol lights up continuously for a short time.
2. If you briefly press the ON/OFF/Mode button **2** when the device is locked, the travel lock symbol will also light up.
3. To deactivate the travel lock, press and hold the ON/OFF/Mode button **2** again for 5 seconds. The travel lock symbol flashes during deactivation. As soon as the travel lock is deactivated, the travel lock symbol goes out. The device starts up.

7. CLEANING AND MAINTENANCE

Clean the device and the brush head **1** after each use. To do so, proceed as follows:

1. Pull the brush head **1** upwards and off.
2. Clean the brush head **1** under running water. Make sure to clean the inside of the brush head **1** as well. **D**
3. Clean the device under running water. However, do not submerge the device in water. Water might enter the device and damage it. **E**
4. Regularly clean the charger **5** with a damp cloth and a mild soap solution. Do not use any abrasive cleaning products. Keep the inductive charger **5** and the USB charging cable with mains adapter away from water.

8. ACCESSORIES AND REPLACEMENT PARTS

To purchase accessories and replacement parts, please visit www.beurer.com or contact the corresponding service address (according to the service address list) in your country. Accessories and replacement parts are also available from retailers.

Designation	Item number and/or order number
Classic Clean brush heads 4-pack	105.21
Premium Whiten brush heads 4-pack	105.31

Premium Plaque Eraser brush heads 4-pack	105.22
Premium Protection Plus brush heads 4-pack	105.23

9. DISPOSAL

For environmental reasons, do not dispose of the device in household waste at the end of its service life. Dispose of the device at a suitable local collection or recycling point in your country. Observe the local regulations for material disposal. Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal. You can obtain the location of collection points for old devices from the local authorities.

The empty, completely flat batteries must be disposed of in specially designated collection boxes, at recycling points or at electronics retailers. You are legally required to dispose of the batteries.

The codes below are printed on batteries containing harmful substances:

Pb = battery contains lead
Cd = battery contains cadmium
Hg = battery contains mercury

Removing the battery

1. To open the electric toothbrush, use the underside of the charger **5**.
2. Place the bottom of the toothbrush on the rotating template on the underside of the inductive charger.
3. Now apply a little pressure and turn anti-clockwise. The end cap on the toothbrush can now be removed.
4. Fold open the bottom of the toothbrush.
5. Pull out the inner housing of the toothbrush.
6. Remove the battery from the inner housing.
7. Dispose of the battery separately in specially designated collection boxes, at recycling points or at electronics retailers.

10. TECHNICAL SPECIFICATIONS

Type	SC 50
Dimensions	Device: 184 x 28 mm
Weight	Device: 161 g (toothbrush + 2 brush heads + charger)
Battery: Capacity Nominal voltage Type designation	800 mAh 3.7 V Lithium-ion battery
Power supply	5.0 V $\overline{\text{---}}$ 1.0 A, 5.0 W

Subject to technical changes

11. GUARANTEE

Further information on the guarantee and guarantee conditions can be found in the guarantee leaflet supplied.

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Beurer GmbH • Söflinger Str. 218
BEURER UK Ltd., Suite 16,