beurer






TB 15 Green Planet

DE	Elektrische Zahnbürste Gebrauchsanweisung.....	4
EN	Electric toothbrush Instructions for use.....	13
FR	Brosse à dents électrique Mode d'emploi	21
ES	Cepillo de dientes eléctrico Manual de instrucciones.....	30
IT	Spazzolino elettrico Istruzioni per l'uso	39
TR	Elektrikli diş fırçası Kullanım kılavuzu	47
RU	Электрическая зубная щетка Инструкция по применению.....	55
PL	Elektryczna szczoteczka do zębów Instrukcja obsługi	64
NL	Elektrische tandenborstel Gebruiksaanwijzing	73
DA	El-tandbørste Betjeningsvejledning.....	82
SV	Eltandborste Bruksanvisning.....	90
NO	Elektrisk tannbørste Bruksanvisning.....	98
FI	Sähköhammasharja Käyttöohje.....	106

CE

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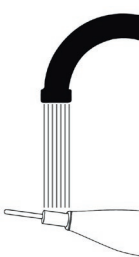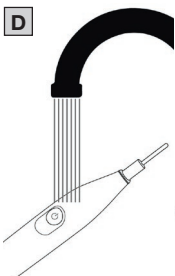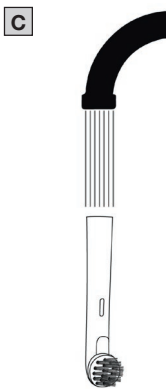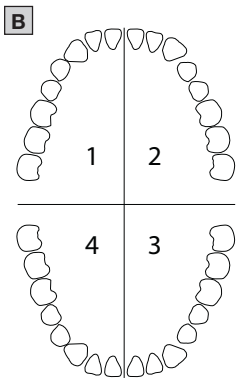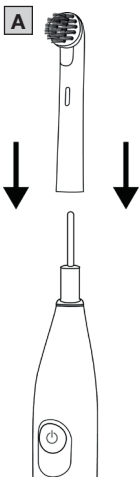ENGLISH

Read these instructions for use carefully. Observe the warnings and safety notes. Keep these instructions for use for future reference. Make the instructions for use ac-

cessible to other users. If the device is passed on, provide the instructions for use to the next user as well.

CONTENTS

1. Warnings and safety notes.....	13
2. Included in delivery	16
3. Signs and symbols	16
4. Intended use.....	18
5. Device description.....	18
6. Use	18
6.1 Preparing for use.....	18
6.2 Starting use.....	18
7. Cleaning and maintenance.....	19
8. Accessories and/or replacement parts	19
9. Disposal.....	19
10. Technical specifications	20
11. Warranty	20

1. WARNINGS AND SAFETY NOTES

⚠ WARNING

- Keep the device and the charging cable away from water or other liquids.
- Do not place the device near a sink or bathtub full of water.

- This device may be used by children over the age of 8 years and by people with reduced physical, sensory or mental capabilities or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely, and are fully aware of the risks of using the device.
- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children unless supervised.
- Do not use the device if it shows signs of damage or does not function properly. In these cases, contact Customer Services.
- If the mains connection cable is broken, you must dispose of the charging cable.
- Only use the USB-C charging cable supplied with your device.

- Never use the device with the mains adapter near a bath, sink, shower or other containers filled with water or other liquids – risk of electric shock.
- The device may only be charged with a mains adapter (not included in delivery) with a safety extra-low voltage output of 5 V and max. 1 A.
- Under no circumstances should you open or repair the device yourself, otherwise proper functionality can no longer be guaranteed. Failure to comply will invalidate the warranty.

Notes on handling batteries

⚠ WARNING

- **Risk of explosion! Risk of fire!** Failure to comply with the points mentioned can result in personal injury, overheating, leakage, venting, breakage, explosion, or fire.
- Always use the correct or supplied charging cable/charger/mains adapter for charging.
- Avoid continuous charging or overcharging. Unplug the charger when charging is complete.
- To charge the device, only use a protection class 2 mains adaptor (), which is approved for use with household appliances and whose output corresponds to the safety extra-low voltages of 5 V and max. 1 A.
- Charge the device under supervision, paying attention to any heat generated, deformation, or release of gases. If in doubt, stop charging.
- Dispose of defective batteries/charging cables/chargers immediately and properly (see chapter on disposal).
- Do not use modified or damaged batteries/charging cables/chargers.
- Always select the correct battery type. Contact Customer Services for information on whether a battery replacement part can be replaced.

- Do not throw batteries into a fire.
- Never forcibly discharge, heat, disassemble, open, crush, deform, encapsulate, modify, or knock batteries.
- Never short-circuit batteries or battery compartment contacts.
- Protect the batteries from direct sunlight, rain, heat, and water.
- Exposure of batteries to an environment with extremely high temperatures or an extremely low air pressure may result in explosion or leakage of flammable liquids and gases.
- Always insert the batteries correctly, taking into account the polarity (+ / -).
- If a battery has leaked, put on protective gloves and clean the battery compartment with a dry cloth.
- If fluid from a battery comes into contact with your skin or eyes, wash the affected areas with water and seek medical assistance.
- **Choking hazard!** Keep batteries out of the reach of children. Seek medical attention immediately if swallowed. Swallowing them may cause severe internal burns and death.
- Never allow children to replace batteries without adult supervision.

⚠ CAUTION

- Store battery replacement parts away from metal objects in a well-ventilated, dry, and cool place.
- Keep batteries clean and dry.
- Keep batteries away from water.

2. INCLUDED IN DELIVERY

Check that the exterior of the cardboard delivery packaging is intact and make sure that all contents are present. Before use, ensure that there is no visible damage to the device or accessories and that all packaging material has been removed. If you have any doubts, do not use the device and contact your retailer or the specified Customer Services address.

- 1 electric toothbrush
- 1 brush head
- 1 USB-C charging cable
- 1 instructions for use

3. SIGNS AND SYMBOLS

The following symbols are used on the device, in these instructions for use, on the packaging and on the type plate for the device:

⚠ DANGER

Indicates an imminent danger. If it is not avoided, it will result in death or serious injury.

⚠ WARNING

Indicates a potentially impending danger. If it is not avoided, there is a risk of death or serious injury.

CAUTION

Indicates a potentially impending danger. If it is not avoided, slight or minor injuries may result.

NOTICE

Indicates a potentially harmful situation. If it is not avoided, the system or something in its vicinity may be damaged.

Product information

Note on important information

Read the instructions

Disposal in accordance with the Waste Electrical and Electronic Equipment EC Directive (WEEE)

Do not dispose of batteries containing harmful substances with household waste

The products demonstrably meet the requirements of the Technical Regulations of the EAEU

CE labelling

This product satisfies the requirements of the applicable European and national directives.

Separate the product and packaging elements and dispose of them in accordance with local regulations.

Marking to identify the packaging material.
A = material abbreviation,
B = material number: 1–7 = plastics,
20–22 = paper and cardboard

Conformity assessed marking for the United Kingdom

Class III device

Importer symbol

Manufacturer

Protection against splashing water

USB / Micro-USB connection

For indoor use only

GRS – standard for independent certification of recycled material

Packaging made from recycled material

4. INTENDED USE

Only use the device to clean teeth on humans. Do not use the device on animals! The device is not suitable for children under 8 years. The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage resulting from improper or careless use.

5. DEVICE DESCRIPTION

The associated drawings are shown on page 3.

1 Brush head

2 ON/OFF/Mode button

Short press: Switches the device on/off

Long press: Changes the mode

3 Light ring

Charge level indicator

Light ring flashes white: Low battery power

Light ring pulsates white: Battery charging

Light ring permanently white: Battery charged

6. USE

6.1 Preparing for use

Charging the electric toothbrush

Fully charge the electric toothbrush before first use.

To do so, proceed as follows:

1. Connect the USB-C charging cable to the electric toothbrush. Plug the USB-C charging cable into a USB port or into a suitable mains adapter and plug it into a socket.
2. During the charging process, the light ring 3 pulsates white. As soon as the electric toothbrush is fully charged, the light ring 3 illuminates permanently in white.

Fitting the brush head

Place the brush head 1 onto the device from above until it clicks into place A. Make sure that the brush head 1 points in the same direction as the ON/OFF/Mode button 2.

We recommend changing the brush head 1 after 3 months. If the bristles of the brush head are worn or frayed, they may need to be replaced more frequently.

6.2 Starting use

1. To switch on the electric toothbrush, briefly press the ON/OFF/Mode 2 button. The electric toothbrush vibrates briefly every 30 seconds to remind you to clean all 4 jaw quadrants for the same amount of time B. The electric toothbrush has an automatic timer and vibrates three times after 2 minutes. The toothbrush switches off automatically after 6 minutes.

2. Select your desired mode by pressing and holding the ON/OFF/Mode button **2**. The following modes are available:

“clean”	For thorough and complete cleaning
“sensitive”	For gentle and gum-friendly cleaning

Differentiation of modes based on rotation speed:
Faster rotation if “clean” is selected
Slower rotation if “sensitive” is selected

3. To cancel a session prematurely, briefly press the ON/OFF/Mode button **2**.
4. Clean the device as described in Chapter “Cleaning and maintenance”.

7. CLEANING AND MAINTENANCE

Clean the device and the brush head **1** after each use. To do so, proceed as follows:

1. Pull the brush head **1** upwards and off.
2. Clean the brush head **1** under running water **C**. Make sure to clean the inside of the brush head **1** as well.
3. Clean the device under running water **D**. Do not submerge the device in water. Water might enter the device and damage it.
4. Keep USB-C charging cable away from water.

8. ACCESSORIES AND/OR REPLACEMENT PARTS

To purchase accessories and/or replacement parts, visit www.beurer.com or contact the corresponding service address (as per the service address list) for your country. Accessories and/or replacement parts are also available from retailers.

Designation	Item number and/or order number
Brush heads Clean Eco Blue 4 pcs	105.27
Brush heads Sensitiv Eco Blue 4 pcs	105.28

9. DISPOSAL

For environmental reasons, do not dispose of the device in household waste at the end of its service life. Dispose of the device at a suitable local collection or recycling point in your country. Observe the local regulations for material disposal. Dispose of the device in accordance with EC Directive Waste Electrical and Electronic Equipment (WEEE). If you have any questions, please contact the local authorities responsible for waste disposal. You can obtain the location of collection points for old devices from the local or municipal authorities, local waste disposal companies or your retailer, for example.

The empty, completely flat batteries must be disposed of in specially designated collection boxes, at recycling points or at electronics retailers. You are legally required to dispose of the batteries.

The codes below are printed on rechargeable batteries containing harmful substances:

- Pb = battery contains lead
- Cd = battery contains cadmium
- Hg = battery contains mercury

Removing the battery

1. Remove the rubber seal in the bottom of the toothbrush with a pointed object.
2. Remove the membrane with a pointed object.
3. The visible screw can now be removed with a suitable screwdriver.
4. Remove the bottom of the toothbrush.
5. Pull out the inner housing of the toothbrush.
6. Remove the battery from the inner housing.
7. Dispose of the battery separately in specially designated collection boxes, at recycling points or at electronics retailers.

10. TECHNICAL SPECIFICATIONS

Type	TB 15
Dimensions	Length: 20.5 cm x diameter 29 mm (width)
Weight	124 g (toothbrush + brush head)
Battery: Capacity Nominal voltage Type designation	600 mAh 3.7 V Lithium-Ion battery
Power supply	5,0 V  1,0 A, 5,0 W

Subject to technical changes

11. WARRANTY

Further information on the warranty and warranty conditions can be found in the warranty leaflet supplied.

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Remarques relatives aux batteries

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Beurer GmbH, Söflinger Str. 218
BEURER UK Ltd., Suite 16,